

# Star<sup>HYGIENE</sup> DENTALEZ

## *BETTER HYGIENE MADE EASY WITH iSTAR™ PROPHY HANDPIECE*

Our iStar prophy features a simple, compact, cordless design that dramatically simplifies the prophylaxis procedure. No more stress from the pull of handpiece tubing. The iStar is ideal for the modern dental operator and features thoughtful design elements that focus on comfort and efficiency.

*Custom disposable barriers for  
added protection for the patient*

*Five speed options from 500-2,500 rpm allow speed  
adjustment based on user preference*

*On/Off fingertip control for intuitive operation*

*Compatible with most disposable prophy angles*

*Part number: 265993*



## PROPHY STAR® 3 AND 3 SW

Constructed with anodized aluminum, the lightweight Propphy Star line helps minimize operator fatigue. The LubeFree design saves time and money with no required maintenance. The 360-degree swivel helps with comfort and operator control.

*Swivel (SW) and fixed back-end configurations*

*Air-driven control up to 5,000 rpm*

*LubeFree design is maintenance free which saves time and money*

*Use with your preferred disposable prophylaxis angle*

*Propphy Star 3 part number: 264422*

*Propphy Star 3 SW part number: 266281* **NEW!**



Specifications	Propphy Star 3	Propphy Star 3 SW
Speed	0-5,000 rpm adjustable	0-5,000 rpm adjustable
Air Pressure	35-40 psi	35-40 psi
Dimensions	0.64" in diameter	0.64" in diameter
Weight	62 grams	79 grams
Length	4.3"	3.71"
Connection	4-line configuration	Swivel

## TITAN® BLIS-SONIC™ SCALER

Use sonic scalers to effectively and gently remove hard calculus deposits and stains. Conveniently connects to existing air-driven tubings and foot control. Ergonomic silicone grips with comfortable positioning reduces hand fatigue.


*Swivel (SW) and fixed (K) back-end configurations*

*Lightweight design*

*Includes five scaling tip designs*

*SW part number: 264563*

*K part number: 264565*

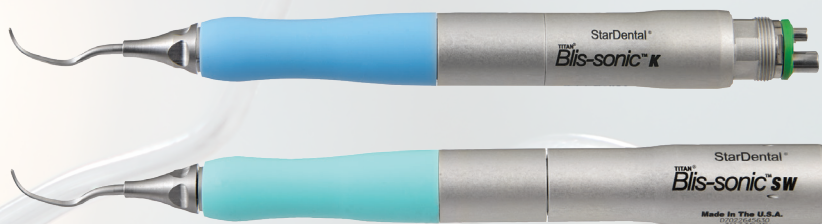

### Star Nomenclature Key

**SW** = Swivel (coupler) back-end

**K** = Fixed back-end

### Blis-sonic SW Specifications

Air Pressure	40 psi
Weight	54.6 grams
Length	6.8 inches (17.3 cm)


### Blis-sonic K Specifications

Air Pressure	40 psi
Weight	50.7 grams
Length	6.7 inches (17.0 cm)



#### Universal

- Use for all surfaces, medium and heavy calculus

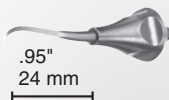

#### Sickle

- Recommended for sensitive patients
- Less powerful than the Universal Tip
- Less suitable for deep root scaling



#### Perio

- Recommended for lingual supragingival scaling of mandibular incisors and deep local pockets
- Provides same energy as Universal Tip
- Also includes left and right perio tips (not shown)